


Fast-Pack pro

Obturation Pack Handpiece



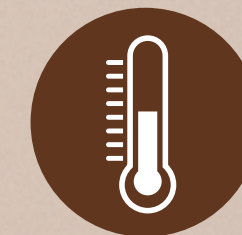

Temperature Range

Wide temperature range of 90°C to 300°C
Five memory programs to customize temperature and cut-off time as per obturation material used or any other intended use of the device



Real-time temperature display

T1 300°C

Internally Heated Tip

Only the working end (4-5mm) of the tip is heated-avoids need for insulating sleeves
The non-working end is insulated and stays cool avoiding soft tissue burn



Class-leading Controlled Heating

Fast heating - Quick heating prevents PDL damage ●

Heating Time 0.2 seconds for 200°C
0.5 seconds for 250°C

Swift cooling - Swift cooling ensure GP does not follow along the heating tip ●

Cooling Time 2 seconds for cooling

Experience Ergonomics

360° start / stop ring button allows holding the device in any angle comfortable for you
Left - right handed operation
Comfortable grip, light in weight



Powerful Battery

Highly efficient 2600mAh Lithium battery
Direct Plug & Play - In case of low battery, use the directly with adaptor device

- Charging Time 2 hours to complete charging
- Operating Time 4 hours of continuous operation at 250°C



Fast-Pack pro

Obturation Pack Handpiece



●
Medium
(Yellow)
50/0.05 - 1 Qty

●
Small
(Black)
40/0.025 - 1 Qty

●
Large
(Blue)
60/0.06 - 1 Qty

